# *Slim de Beauté*

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Instruction Manual

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**CL-BT-10**

# *Slim de Beauté*

## Instruction Manual

Please read this manual beforehand

- Thank you for your purchase.
- Please read this instruction manual beforehand, and use it appropriately. Please make sure to read "Safety Notice" before using it. Then use the device safely.
- Please keep this instruction manual with great care.
- This device is for cosmetic purpose. Medical use and sales are prohibited in a law.
- It expectes to be applied to the healthy skins. We are not in a position to take any responsibilities for any troubles to the use to the part treated with a doctor.

# Safety Notice

- Please read the “Safety Notice” carefully beforehand, then use the device appropriately. Also, please keep the manual with great care in case you need it.

To prevent human hazards and damages of possessions beforehand, please read these instructions

## Safety Notice

Belows are explained the classifications of hazards and damages with images after you ignored these instructions and mistakenly operated. Please understand these before you read the manual.

### **Danger**

It threatens death or sever injuries

### **Warning**

It could be in danger of death or sever injuries

### **Attention**

It could lead to getting injuries or damaging possessions

## ■ Means of images below

Marks explain what to be prohibited

Mark explains what to be obligated

Mark explains what to be taken care for

## Warning

	Please read this manual beforehand.
	Please stop using it and consult with a doctor if any abnormalities happen.
	Please turn off and contact to the purchased shop or origin of the seller if abnormal sound, smells, smokes, or allotrio were found.
	This product expects for body use only . Never apply it on face and nearby heart .
	Stimulation excites muscles. Make sure no dangerous materials nearby.
	Stimulation could concentrate on one part after the viscosity of the gel pad is losen. Please follow the target indication of use.
Please discard the battery properly. Never throw into a fire.	
	Never operate it by wet hands. It threatens shock, breakdown, fire disaster, and burn.
	Never soak into water and clean it with water .
	Never put allotrio inside.
	Never break, repair, and remodel it.

## Attention

	This device is used for household. Never apply it for institutional use.
	Never apply it to small kids.
	Never use while driving and sleeping.
	Please avoid metals like rings, watches, glasses and other accessories from your body while using it.
	Please avoid to drop the device when using it. It could cause o f break.
	Avoid water. Never use and keep it in a bath, hot, humid , and dusty places.

# Other Attention


Never use it if you are in these conditions.

## Warning

- 1, Heart disease
  - 2, Life support medical electric device like an artificial heart -lung, mounting type of electric device like an electrocardiogram, and intracorporeal implanting type of electric device like a pacemaker
  - 3, Abnormal blood pressure
  - 4, With a high fever
  - 5, Parts with a fever, sudden illness, skin diseases, keloid constitutions, high blood pressure, during menstruation, expected mothers, after accidents, stomach with myoma of the uterus, extraordinary fatigue, conditions when resting is necessary
  - 6, Body parts with metals, plastics, or silicon inside
  - 7, Never use it when you feel abnormal
  - 8, Never use it to the sensitive areas like near heart, scratches, moles, vein, inside the mouth, and the genital areas
  - 9, Atopic dermatitis or sensitive skins with allergy
  - 10, Acute diseases, transmissible diseases, and infectious diseases
  - 11, Being anemic, diabetes, hemophilia, and tuberculous disease
  - 12, Long term use of steroid-based hormone drug or telangiectasis by the liver function failure
  - 13, Under drinking
  - 14, Under medications or taking medicines
  - 15, Abnormal skins like an eruption by cosmetics, inflammation, sunburn, itchy skins, scratches, possibility of bleeding, and perceptual dysfunction
  - 16, Abnormal, inflammation, or scratch in your breasts
  - 17, Please consult with a doctor if you get plastic surgeries or breast enlargement surgeries
  - 18, Under 15
  - 19, Maligrancy
  - 20, Taking or injecting painkillers, hypoglycemic drugs, or hypotension drugs
  - 21, Please avoid to apply it to stomach just after meals
- ※ Please consult with a doctor if you regularly go to hospitals
- ※ Please stop using it and consult with a doctor if you are worried about using it or had abnormal while using it.

# Product Specification


Body × 1

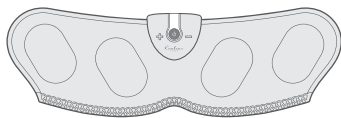Safekeeping Tray × 2

Gel Pad × 2

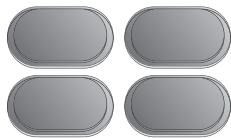※Sale for spare  
Spare gel pad 4 pieces 2,800yen +tax

Alkali AAA battery × y 2

【For testing】

Standard name: Slim de Beaute

Classification: Household EMS exercise apparatus

Model number: CL-BT-10

Body size: W427mm × D140mm × H22mm

Body weight: 122g (body only)

Battery type: Alkali AAA battery × 2

Material

【Body】Silicon, ABS plastic, conductive silver ink

【Gel Pad】Hydrogel

【Safekeeping Tray】PET


Components

Body, Gel Pad × 2, Safekeeping Tray × 2,

Alkali AAA battery × 2 (For Testing),

Instruction Manual (Guarantee Attached)

# Get Started
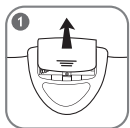
## ■ Set the battery

① Remove the battery cover


② Set two alkali AAA batteries as instructed image. Close the battery cover afterwards.

※ When there is a gap between the main body and battery cover, sweat invades the battery cover inside and causes the trouble. Please fasten the battery cover.

## ■ Attach gel pads

The viscosity gets loosen and consumes faster than expected in case you apply inaccurately

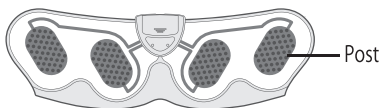The blue film	Sticks to the posts in device
The clear film	Sticks to your skins

- ① Remove the blue film first.
- ② Attach all the gel pads to each posts horizontally so that air does not enter if possible.
- ③ Stick them well with ball of a finger to avoid air inside from the clear film.
- ④ Remove the clear film and get ready.

# How to use

## 1 Attach the main body

Attach the main device to the parts you like to get treatment referring to "Main device mounting location".

※ Not turned on electricity when the gel pads do not stick fast to your skin.


## 2 Turn on the device

Treatment starts after pressing the power button for 2 seconds, and sound, turn on then.

※ It starts from level 1 when you turn on.

## 3 Adjust the level

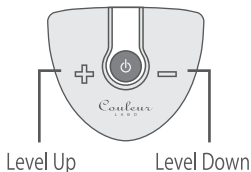Short press "+" and select the level. Then a sound sounds.

※ The output level does not increase if you press the button too long.

Short press "-" to decrease the level. Select your comfortable level and start treatment. It turns off automatically after 20 minutes.

Press the power button for 2 seconds to turn off in mid.

15 levels are available

## 4 Finish

It automatically turns off after 20 minutes. Please remove the device to the vertical direction carefully. Then wipe off sweat and moisture, and storage it attached safekeeping film.

Safekeeping film


Please do not bend it and pull it when you storage it.


## Mode auto-switch


3 modes change in return for 20 minutes to get effective exercise result.

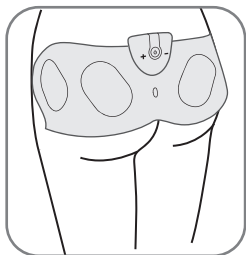
### Training Pulse
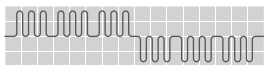
### Smoothing Pulse
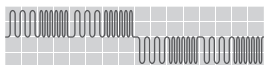
### Firming Pulse
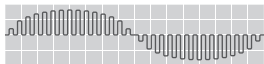
## Precautions for use

- Please stop using it and consult with a doctor if you feel abnormal.
- Even if you use it with the stronger level and for a long time, it does not promise to rise the effects. It could cause of the muscle fatigue or skin troubles.
- Please avoid metals like necklaces or other accessories from the gel pads.
- Please do not use body cream, lotion, and oil during the treatment.
- It is easier to stick fast to your skins if body hair is removed.
- Please do not pull the device forcedly when you use it.
- Please avoid to press buttons too strong.
- Please do not sit down while you are applying the device.

# Main device mounting location

## Hip

Apply it at the center of your coccyx.

## Waistline

Apply it at the center of your backbone

## Stomach

Apply it at the center of your bellybutton.

※ It works only if it is applied to your skins properly.  
The electric stimulation does not happen if the device is not applied to your skins.

## Recommended treatment

### 4 tips for effective results

- ① Take a break once or twice a week.
- ② Record the maximum of power level daily, increase the level gradually.
- ③ The stronger the output level could lead to a better result, but please avoid the excessive stimulation so as to be unpleasant.
- ④ Continue the treatment until it automatically turns off for 20 minutes.


## ⚠ Must Follow

※ Take a 60 minutes break in the same part of your body , and apply two times a day in maximum.



# Daily maintenance & Attentions

## ⚠ Attention

- Please make sure to turn it off when you clean it.  
You could receive shock if it is cleaned with wet towel or tissues.

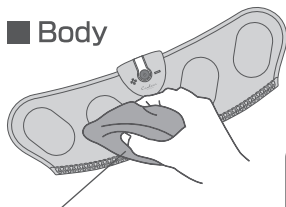Dry towel or towel with some soapy water

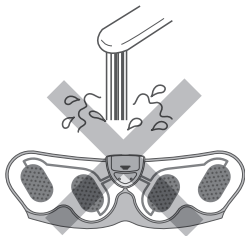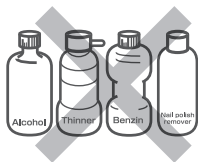- Never clean it with water.
- Please do not use thinner, benzine, and alcohol.

It could cause of break, cracking the parts, and changing color.

Please wipe off with dry towel or towel with some soapy water.

## ■ Gel Pad

- Never attach sticky sides to each other.
- Please store the gel pads attached to the sticky side in the attached tray.
- Please clean the gel pad with a drop of water to a finger-tip.
- Please keep it in an airy place without direct sunlight.

## ■ Battery

- It is indicated about 20 times available for 20 minutes treatment a day.
- ※ It depends on operating time and the usage count.
- Please make sure to turn off when changing the battery.
- Please set CR-2032 round type of lithium battery.

Make sure to set the battery properly.

- Please keep it in an airy place without direct sunlight.
- To keep it last longer, please remove the battery from the device when you do not use it. It could get corrosion on its face if it is kept inside the device long time.
- Please discard the battery properly. Never throw it into a fire.

## Consumable

- ※ The gel pads are consumable. It expects to be used 20 times. Please change the gel pads if the viscosity gets losen. Please contact to the purchased shop or the origin of the seller if you have any questions.

## ■ Spare Gel Pads 8 pieces

Model number: CL-BT-10PAD 3,600yen + tax

# Storage

## ■ Please keep the device in these places

Flat & Stable place	Airy place	Place of the temperature in the image
		

## ■ Never keep the device in these places.

This device is made of accurate machines and electric parts.  
It could cause of malfunction if you keep it in following areas .  
Please avoid the following areas.

With direct sunlight 	Lots of dust 	Thermal change 
Humidity change 	Place with the fire 	Humid 
Place with the volatile 	Near the air-conditioner 	Place with the vibration 

## Attention

- Please do not put it near small babies and persons who cannot express feelings in words.  
It could cause of accidents and deconditioning.

# Troubles Q&A

	Procedure
Stimulation to be hurt	<ul style="list-style-type: none"> <li>• The gel pads are consumed/not stick to skins fast Please change the gel pads Please make sure to stick to skins firmly</li> <li>• Muscles get tired Must follow the instruction of use, twice a day in the same part of your body Please take 60 minutes break when to apply to the same part of your body.</li> </ul>
The gel pads stick to skins too tight	<ul style="list-style-type: none"> <li>• After applying the gel pads in the posts, please stick on the device well to let air remove.</li> </ul>
Change color to brown or like black around the posts (silver area)	<ul style="list-style-type: none"> <li>• Sweat and other moisture cause this. It does not affect its function and usage.</li> </ul>
<p>The output level does not increase</p> <p>The error sound does not stop</p> <p>Does not work after turn it on</p> <p>Does not feel energization</p> <p>Does not feel treatment</p>	<ul style="list-style-type: none"> <li>• Operation error Make sure not to long-press "+", or to press "+" firmly.</li> <li>• The batteries are consumed Please change new batteries and check it again.</li> <li>• Warning sound of non-electric state sounds Make sure the gels are stuck to skins firmly. Make sure to remove the clear film.</li> </ul>
Reddening skins where treatment was applied	<ul style="list-style-type: none"> <li>• Stop treatment immediately. Avoid using it until your skin returns normal. Please consult with a doctor if it does not heal</li> </ul>
Loose the viscosity and not stick to skins	<ul style="list-style-type: none"> <li>• Clean the gel pads with a drop of water on your finger-tip.</li> <li>• Must storage the device stuck in the safekeeping sheets . The viscosity extremely gets lost if it is stck in other areas.</li> <li>• Change the gel pads</li> </ul>

# Guarantee & After-service

## ● Guarantee

The guarantee is attached on the end of this manual.

Please make sure that you fill in the information about purchased shop and the date of purchase, and keep it after you read the manual carefully.

The warranty period is 1 year after your purchase.

This guarantee covers only the main body, but not the attached tray, the gel pads, the lithium battery, and other consumable items.

## ● Before asking for repair

Please read this manual and kindly check your device again before you ask for repair. If you have any problems, please contact to the purchased shop.

## ● Breakdown

Please contact to the purchased shop and ask for repair.

## ● Repair after its warranty period is over

Please contact to the purchased shop. If it is available to repair your broken device, the repairing cost will be charged.

## ● The minimum possession period of covered item

The repairing parts of this device is expected to be retained at least 5 years. Therefore, the repairable term of this device is 5 years in principle.

It could be repaired after the period is over, please contact to the purchased shop then.

## ● Important note in case you need a repair

In case you ask the purchased shop or us, repairing cost less than 10,000 yen will be repaired as usual. The cost over 10,000 yen will not be repaired until the purchased shop or the customer gets contacted and agrees

**1. We will repair the broken device for free subject to the device is used under the instructions written in the manual and still in the warranty period. If you need the free repair service, please ask to the purchased shop with the device and this guarantee.**

**2. We will charge for repairing cost under these conditions.**

- ① Break or damage because you use it for other purposes except for written purposes in the manual.
- ② Break or damage because of your unreasonable repairing or remodeling.
- ③ Break or damage because of storage places, transportation, or falling after your purchase.
- ④ Break or damage because it is used for other purposes like an industrial use other than household use.
- ⑤ In case of no indication in the guarantee.
- ⑥ In case of no date of purchase, name, and the purchased shop in the guarantee.
- ⑦ Break or damage by loading on vessels or vehicles.
- ⑧ Fire disasters, earthquakes, flood damages, lightning damages, and other unexpected phenomenon, abnormal voltage, and unspecified power.
- ⑨ Break or damage because you had the device repaired by other than us.

**3. This guarantee is available only in Japan.**

**4. This guarantee is not reissued. Please keep it carefully .**

- This guarantee promises free repairment subject to the period and conditions written in this guarantee.  
Therefore, this guarantee does not declare the customer's judicial right. Please contact to the purchased shop if you have any questions about the repair service after the warranty term is over .
- Please find the after-service in the manual about the repairment after its warranty term is over or the retain period of repair parts

In case you discard the device, please discard as non-burnable refuse as follows the rules of each local governments.

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